




STAIN

POLYFIX

CERAMIX

CLASS CERAMIC

Dental Full-Spectrum Color Solutions

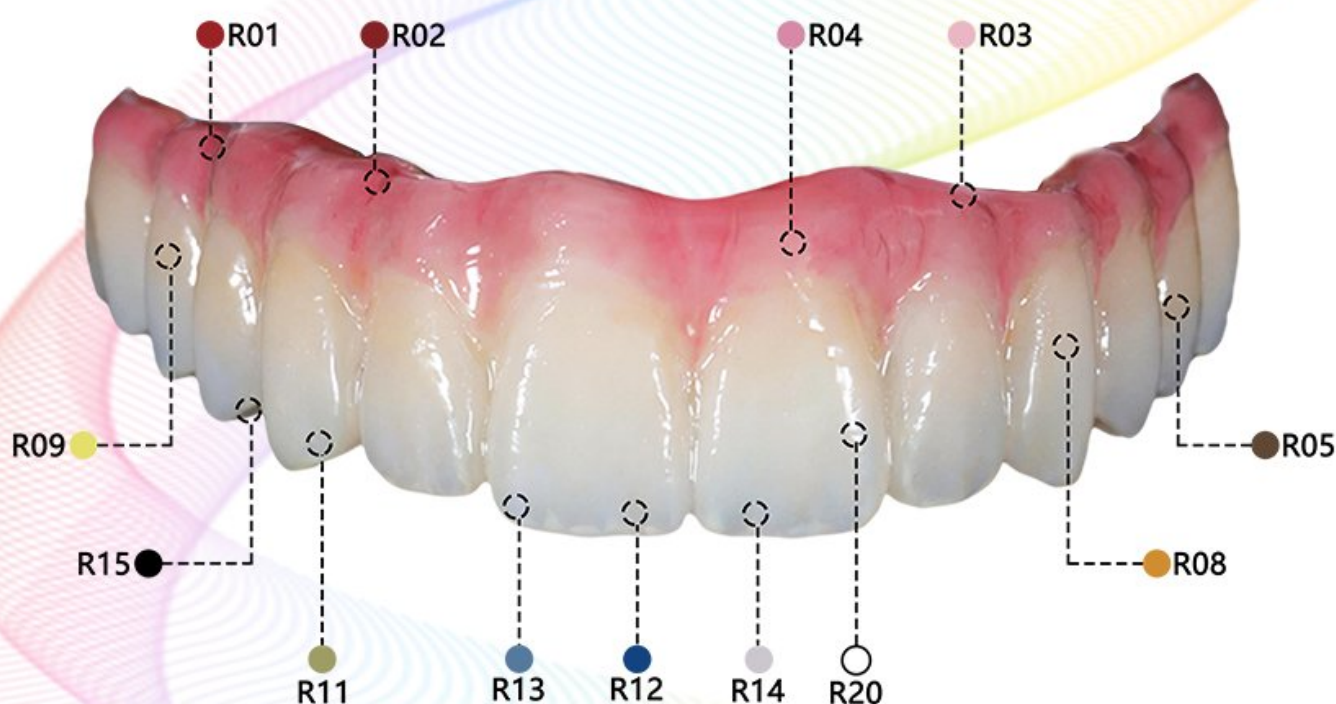

<https://www.zonmedtech.com>

COMPANY PROFILE

5000+


Cooperative Labs & Clinics

2019

Establishment

95%

Customer Retention Rate

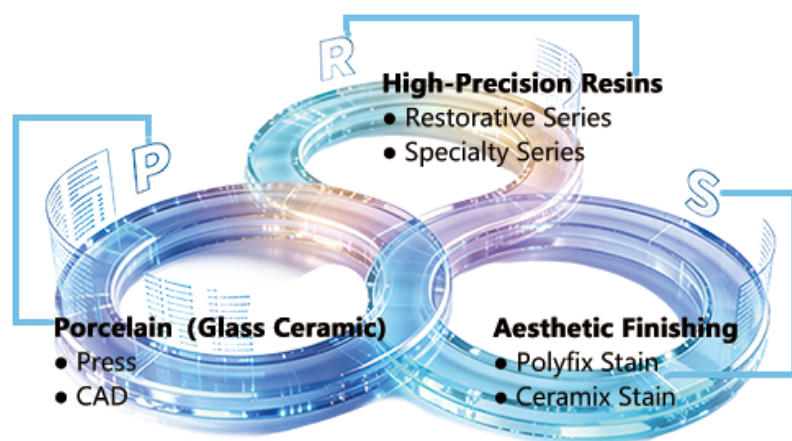

PRS™ ECOSYSTEM

Closed Loop: Materials to Aesthetic Restoration

Founded in 2019, ZonMed is a leading high-tech enterprise dedicated to the research, development, and production of the full PRS Ecosystem.

ZonMed PRS™ provides a seamless, closed-loop digital workflow from initial scan to final restoration. With our own advanced manufacturing facility, we ensure strict quality control and continuous innovation.

We serve the global market directly, exporting our reliable and aesthetically superior solutions to dental professionals and partners in over 60 countries and regions. ZonMed is committed to being your trusted partner for high-performance restorative materials that deliver consistent results worldwide.



POLYFIX STAIN KIT

- Kit A (20 btl)..... P01-02
- Kit B (7 btl)..... P03-04
- Kit A (20 syr) PRO..... P05
- Glaze P06



CERAMIX STAIN KIT

- Kit C (20 btl)..... P07-08
- Kit C (20 syr) PRO..... P09



GLASS CERAMIC

- Press..... P10-11
- Cad..... P12-13
- Sintering Paste..... P14
- Electric Airbrush P15
- UV Curing Light P16




01

POLYFIX STAIN KIT(A)



No.:PRS-S-R(A)

This set achieves natural appearance restoration and personalized aesthetic customization of restorations such as removable denture bases, temporary bridges, and crowns through precise color matching and creative color shaping.

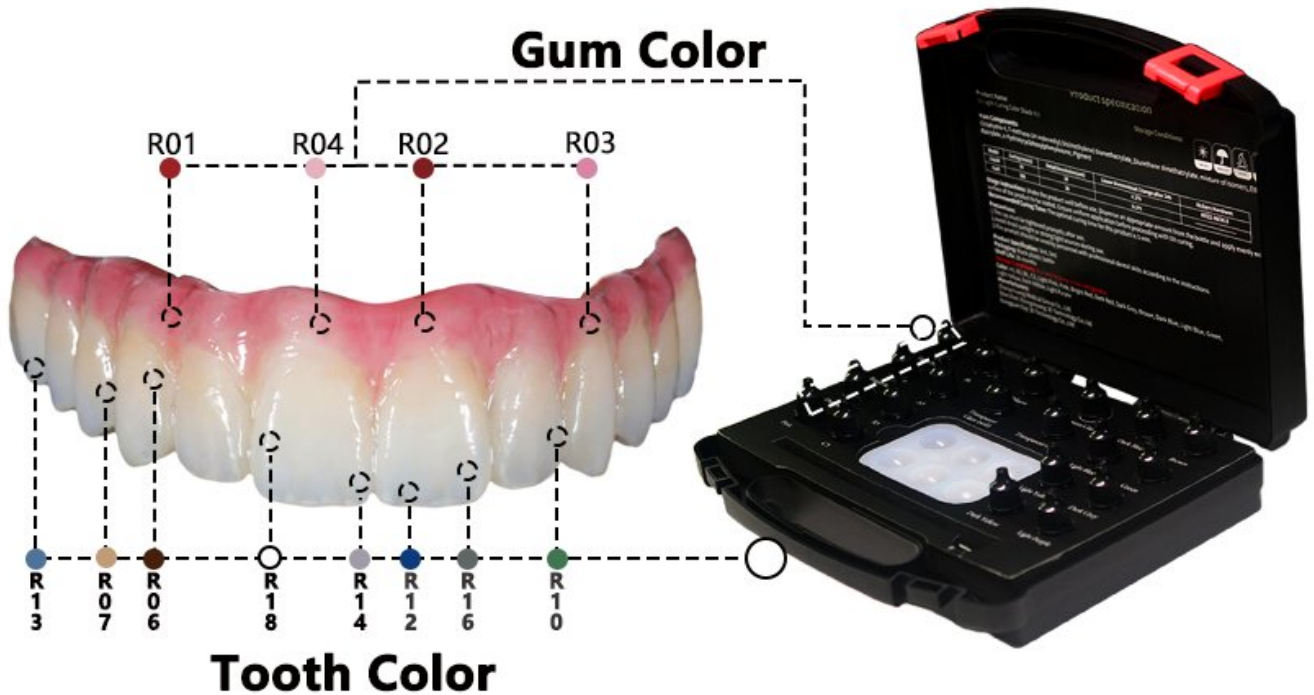

PRODUCT FEATURES

- High Color Stability.
- Fast Curing Speed($\leq 30s$).
- Multi-shade Combination.
- Intelligent Color Matching Technology Nanomaterials with Certain Antibacterial, Wear-resistant, and Corrosion-resistant Properties.

PRODUCT DISPLAY



03

POLYFIX STAIN KIT(B)



No.:PRS-S-R(B)


This set achieves natural appearance restoration and personalized aesthetic customization of restorations such as removable denture bases, temporary bridges, and crowns through precise color matching and creative color shaping.



PRODUCT FEATURES

- Antibacterial and Durable.
- Convenient Post-maintenance.
- Wear-resistant and Scratch-resistant Easy Operation.
- Make the Product Transparent and Glossy.
- Fast Curing Speed($\leq 30s$).

04

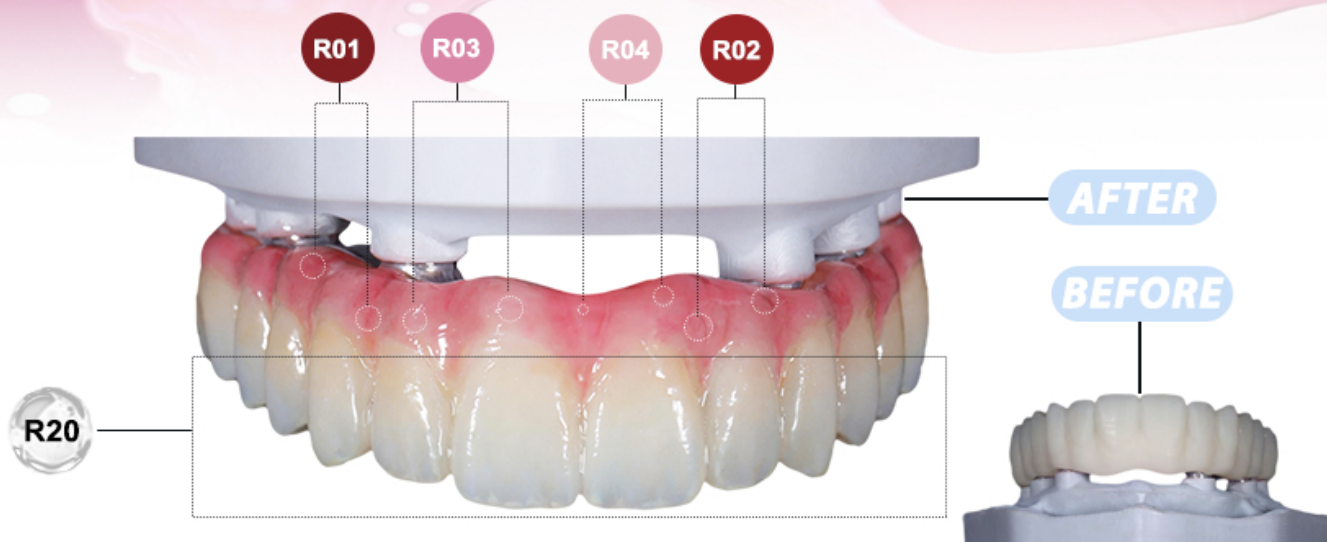

PRODUCT DISPLAY



Application : Gum, Full Mouth Movable Base, PMMA Block, Peek, and other Polymers.

05

No.:PRS-S-R (PRO)

PRODUCT DISPLAY



Bright Red

SR01



Dark Red

SR02



Pink

SR03



Light Pink

SR04



Sepia

SR05



Brown

SR06



Khaki

SR07



Orange

SR08



Yellow

SR09



Green

SR10



Dark Green

SR11



Dark Blue

SR12



Light Blue

SR13



Lilac

SR14



Black

SR15



Gray

SR16



Beige

SR17



Diluent

SR18

TRANSPARENT GLAZE

SR20 & SR21(SOFT)

PRODUCT FEATURES

- Long-lasting Shine.
- Fast Curing Speed ($\leq 30s$).
- Excellent Adhesive Durability.
- Nanofilled Light Cured Protective Coating.



Color Transparent
Net 5ml / 15ml

APPLICABLE PRODUCTS(SR21)

- Resin: Flexible Base, Smile Veneer, Night Guard Pad, Clear Aligner
- Other: PMMA

DETAIL INFORMATION

Transparent Glaze excels at filling tiny surface voids, enhancing smoothness. This glaze not only improves wear resistance and aesthetics but also extends the lifespan. By shielding from the oral environment, it ensures these dental items maintain optimal performance and appearance, serving users better.


AFTER**BEFORE**

07

NO.:PRS-S-G

CERAMIX STAIN KIT

Full series of 20 standard colors, unlimited free collocation

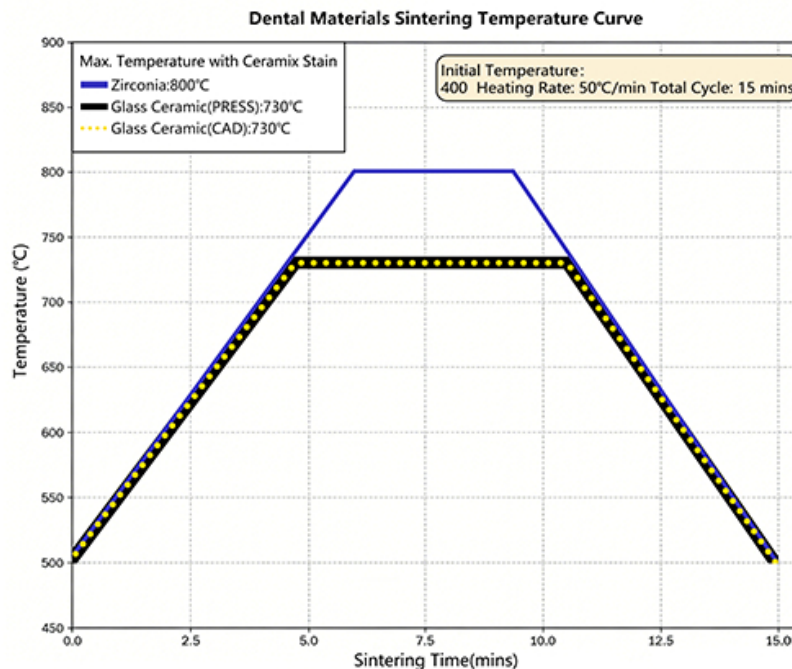


CORE ADVANTAGES

Universal Compatible with all ceramic materials and multiple aesthetic staining techniques.

Flexible Various combinations mix to create any shade.

Easy to Use Innovative paste formula for hassle-free application.

Realistic Stable high-gloss finish (no whitening/graying). Visible natural mimicry.

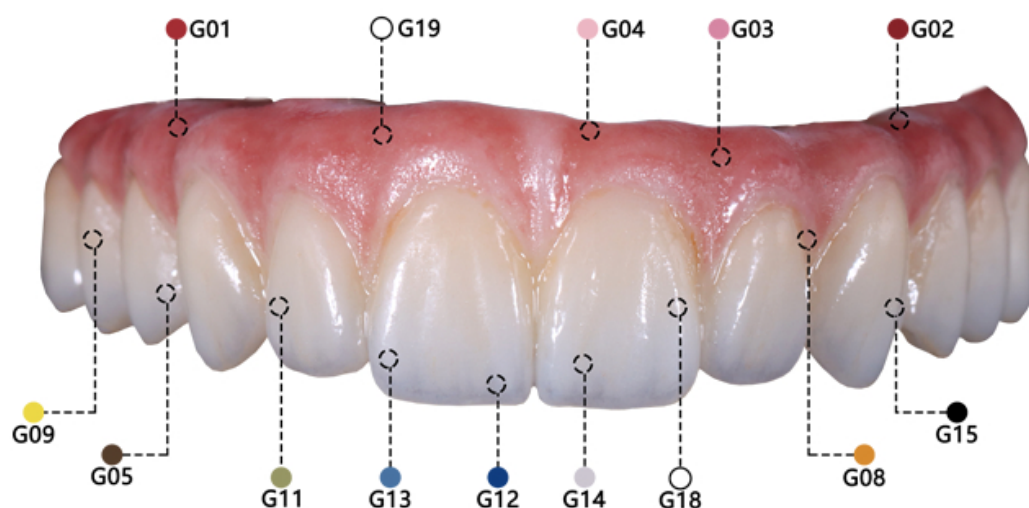

PRODUCT DISPLAY



No.:PRS-S-G (PRO)

PRODUCT DISPLAY



Bright Red Dark Red Pink Light Pink Sepia Bronze Khaki Orange Yellow Green

G01

G02

G03

G04

G05

G06

G07

G08

G09

G10



Dark Green Dark Blue Light Blue Lilac Black Gray Beige White Glaze (FP) Glaze

G11

G12

G13

G14

G15

G16

G17

G18

G19

G20

GLASS CERAMIC(PRESS)

PRODUCT HIGHLIGHTS

- Natural Color Shade, Excellent Aesthetic Effect.
- Minimal Reaction Layer, Easier Clinical Operation.
- Reliable Strength, no Chipping Risk.
- Ideal For Ultra-Thin Veneers, Superior Acid and Stain Resistance.



Heat-pressed ceramic is the ideal material for ultra-thin dental veneers, with core advantages of natural color, acid resistance, no chipping risk, easy operation and excellent aesthetic restoration effect.




Clinical Case 01




Clinical Case 02

COLOR CODE GUIDE

LT

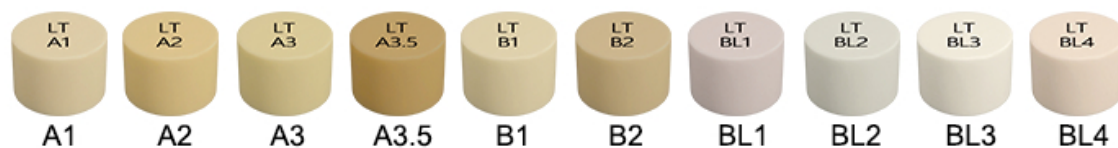

HT

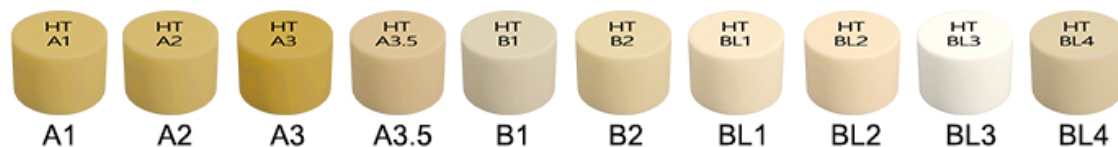

TECHNICAL DATA

Property	Density	Thermal Expansion Coefficient	Biaxial Flexural Strength
Specification	2.4–2.6	$(10 \pm 0.5) \times 10^{-6}$	≥ 420
Unit	g/cm ³	K ⁻¹	MPa

Property	Sintering Temperature	Chemical Solubility	Radioactivity
Specification	920±10	< 100	≤ 1.0
Unit	°C	μg/cm ²	Bq/g

GLASS CERAMIC(CAD)

A high-performance lithium disilicate glass ceramic with superior strength, natural aesthetics and precise processability for dental restorations.

PRODUCT HIGHLIGHTS

- High Strength & Toughness.
- Natural Aesthetic Effect.
- Precise CAD Processing.
- Wide Clinical Applicability.
- Excellent Bonding Performance.



TECHNICAL DATA

Property	Density	Density(After Sintering)	Biaxial Flexural Strength
Specification	2.4~2.6	2.4~2.7	≥480
Unit	g/cm ³	g/cm ³	MPa

Property	Thermal Expansion Coefficient	Sintering Temperature	Chemical Corrosion Resistance	Radioactivity
Specification	$(10 \pm 0.5) \times 10^{-6}$	820±10	<100	≤1.0
Unit	K ⁻¹	℃	μg/cm ²	Bq/g

COLOR CODE GUIDE

H

									
A1	A2	A3	A3.5	B1	B2	B3	B4	C1	C2
									
C3	C4	D2	D3	D4	BL1	BL2	BL3	BL4	

L

									
A1	A2	A3	A3.5	A4	B1	B2	B3	B4	C1
									
C2	C3	C4	D2	D3	D4	BL1	BL2	BL3	BL4

M

									
A1	A2	A3	BL1	BL2	BL3	MOX0	HOX0	MOX1	

SINTERING PASTE

Suitable for veneers, inlays, onlays, partial crowns, single crowns, bridges and other restorations that are inconvenient to place directly on the sintering tray or sintering pin during porcelain baking or glazing processes.



- **Stabilize Fixing**
Secures hard-to-hold restorations during porcelain furnace sintering.
- **Prevent Discoloration**
Restorations stay colorfast post-firing.
- **Easy Removal**
Peels off easily post- sintering, no sandblasting, zero residue.
- **Save Drying Time**
No extra pre-drying in the porcelain furnace.



15

UV CURING LIGHT

Clear In Daylight, Professional Black Filter



BRAND: **ZONMED**

WEIGHT: **35g**

POWER: **40W**

CHARGING: **USB**

MATERIAL: **High-strength Aluminum Alloy**

Note: All data on this page is factory-tested. Subject to actual product in case of deviation or upgrade.

ELECTRIC AIRBRUSH

Designed for spraying glaze on dental models. Easy to operate with fine atomization, it enhances model gloss and extends service life, fitting various dental teaching, display and restoration models.

Ultra-fine Nozzle
for Detailing



Wireless Use
Anywhere
Anytime

Powerful &
Consistent
Spray



For Model
Restoration &
Brightening



Even Paint
Smooth Coverage



Non-stop
Painting &
Touch-ups



Zhongchuang Medical Group Co., Ltd.
Shenzhen Zhongchuang 3D Technology Co, Ltd.



Whats App



WeChat



✉ info@zonmedtech.com

🛒 <https://www.zonmedtech.com>



B10, RM 2703 NEWTREND CENTRE, 704PRINCEEDWARD ROAD EAST, SAN POKONG, KLN, HK

Room 204, Administration Building, Jia Tai Pharmaceutical,
No.1 Industrial Zone, Tangwei Community, Fuhai Street, Bao'an District, Shenzhen City